

### Safety. Repeatability. American Made.



# **Screen Gaskets**





### **Screen Gaskets**

Fluid Path's sanitary screen gaskets are manufactured using FDA and USP Class VI compliant elastomers and fluoroelastomers, with stainless steel screens or mesh. Typically used to filter particulates from liquids, FPP's screen gaskets are easy to install and provide reliable, durable performance. Tri-Clamp configurations are standard, however screen gaskets can be customized on request.

### **Specifications**

Gasket Materials Available: FKM, Buna, EPDM, PTFE, Platinum Silicone, Metal Detectable Elastomers Sizes Available: 1/2" - 6" for Tri-Clamp connections. Other styles are available on request. Screen Material Available: 304 or 316L Stainless Steel Mesh Sizes Available: 10-100 Temperature Range: -50°F to 500°F (Varies by Gasket Material selected) Durometer Hardness: Shore A: 30A - 70A Colors Available: Black, White, Transparent

#### Advantages of some gasket materials:

Buna: Excellent resistance to fats, oils, and petroleum-based fluids FKM/Viton: Acid resistant, good chemical stability, high temp-resistant EPDM: Excellent cell culturing & bioprocess compatibility. Excellent steam resistance PTFE: Highly inert to most chemicals, food, beverage and bioprocess media. Wide temperature range Platinum Silicone: Highest Purity and high temperature resistance for demanding bioprocess applications Flouro-Steel: Highest temperature resistance, does not cold flow, chemically inert.






## Why Fluid Path Products?

Fluid Path Products specializes in quality manufacturing and innovation. To ensure particulates are effectively filtered from process media, FPP's stainless steel screen gaskets are made using only the highest purity, top-quality materials. With a comprehensive range of connection sizes, materials, and screen mesh sizes, fluid filtration is simple, efficient, and adaptable to various applications within the food and biopharmaceutical process industries. The Fluid Path team of in-house engineers and experts are constantly adapting to new situations and implementing designs for the best possible solutions.

All of our gaskets are manufactured in New Jersey and Texas.

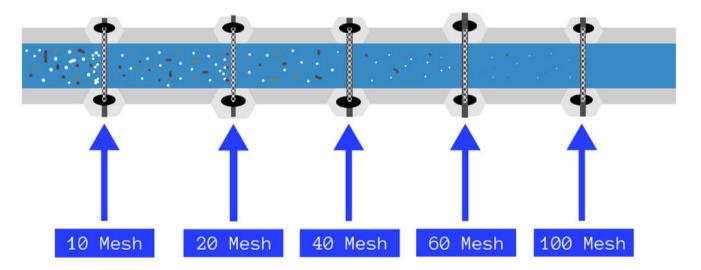

Fluid Path Products customizes! Our in-house tooling design can take your ideas from concept to finished product, typically, within two to three weeks at a cheaper cost than our competitors.